

LiNA MaxFlow™ & LiNA MaxFlow™ Endohook

Powerful suction and irrigation system
with electrosurgery applications



LiNA MaxFlow™ & LiNA MaxFlow™ Endohook

Designed for suction and irrigation.

LiNA MaxFlow™

Suction probe with relieving eyes designed for procedures where delicate suction is required.

LiNA MaxFlow™ Endohook

Monopolar endohook with safety function for retraction of electrode tip into cannula.

LiNA MaxFlow™

Variety of lengths, electrode tips and surgical applications.

Suction probe available in 5 mm or 10 mm, and 5 mm extended length (370 mm).

Choice of solely suction and irrigation function or supplement with electro-surgical monopolar accessory provided by LiNA MaxFlow™ Endohook.

LiNA MaxFlow™ Endohook

L-shaped tip on monopolar endohook.



Safety on/ready
for suction
and irrigation



Safety off/ready
for electro-surgery


Click-on locking mechanism



LiNA MaxFlow™
Endohook



Relieving eyes

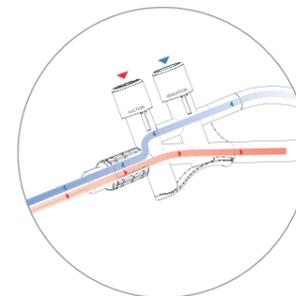

Separate channel
design

Beneficial features to meet surgical needs.


Unique separate channel design in handle designed to reduce risk of blood clots or debris blockage.

True 9 mm direct flow through whole suction system when using a 10 mm probe during suction and irrigation procedure.

No need for frequent change of instruments during surgical procedure as suction, irrigation, dissection, coagulation and cutting is all in one instrument!



LiNA MaxFlow™ & LiNA MaxFlow™ Endohook

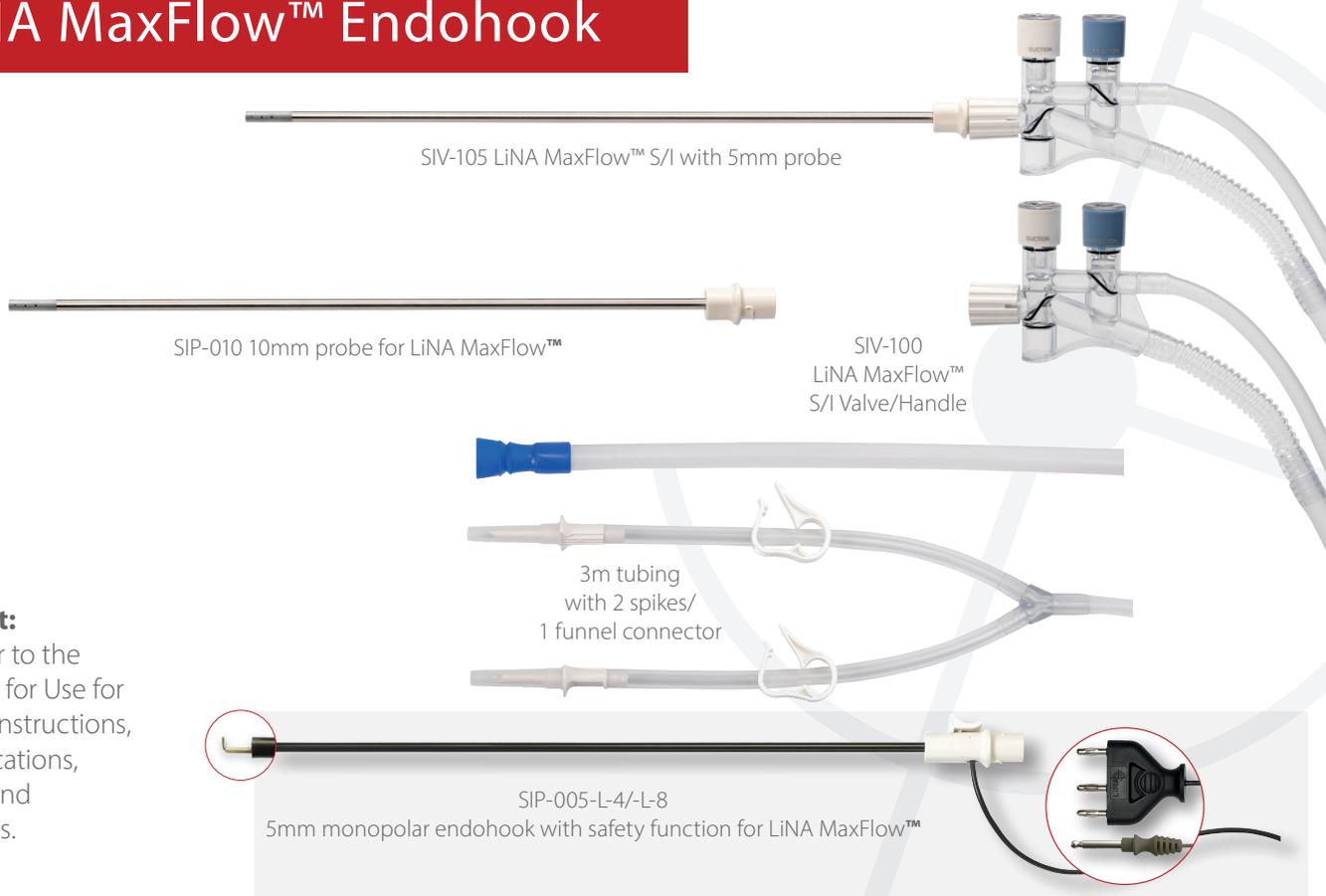

Important:

Please refer to the Instruction for Use for complete instructions, contraindications, warnings and precautions.

Ordering information:

LiNA MaxFlow™ and LiNA MaxFlow™ Endohook



Product code	Description	UOM
SIV-100	LiNA MaxFlow™ (valve only) – 3m tubing and 2 spikes plus connector. Sterile, single use.	Box of 10 ea
SIV-105	LiNA MaxFlow™ suction/irrigation set – with 5mm probe, 3m tubing and 2 spikes plus connector. Sterile, single use.	Box of 10 ea
SIP-005-370	5 mm probe, extended length (370 mm). Sterile, single use.	Box of 10 ea
SIP-010	10 mm probe. Sterile, single use.	Box of 10 ea
SIP-005-L-4	LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 3-pin connector, L-shaped tip electrode. Sterile, single use.	Box of 10 ea
SIP-005-L-8	LiNA MaxFlow™ Endohook – Endoscopic monopolar hook for LiNA MaxFlow™, 3m cable, 8 mm connector, L-shaped tip electrode. Sterile, single use.	Box of 10 ea

EU MDR 2017/745 devices class:

- IIa LiNA MaxFlow,
- IIb LiNA MaxFlow Endohook.

Name and number of the Notified Body: TUV Nord Polska Sp. z o.o. # 2274.

©2025 LiNA Medical AG. All rights reserved.

 LiNA Medical AG, Platz 3, 6039 Root D4 | Switzerland
Phone: +41 78 617 56 56, info@lina-medical.com | www.linamed.com

 LiNA Medical Polska Sp. z o.o., Rolna 8A, Sady, 62-080 Tarnowo Podgórze | Poland
Phone: +48 61 222 21 21, info.poland@lina-medical.com | www.lina-medical.pl

EN ISO 13485:2016
Quality Management System

